


# -30°C Biomedical Freezer

DW-30L698

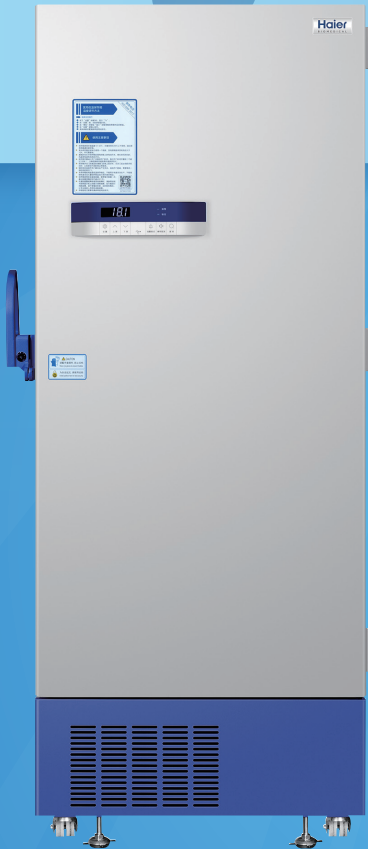

## Scope of Application:

Ideal for utilization in blood banks, hospitals, disease control centers, research institutes, as well as the electronic and chemical industries. Specifically designed for the cryopreservation of plasma, biologics, and other related products.



### Qingdao Haier Biomedical Co.,Ltd.

No.280 Feng Yuan Road, High-tech Zone,  
Qingdao, 266109, P.R. China  
E-mail: [inquiry@haierbiomedical.com](mailto:inquiry@haierbiomedical.com)  
Website: [www.haiermedical.com](http://www.haiermedical.com)



Haier Biomedical  
International



Haier Biomedical  
International



@haiermedicalint



Haier Biomedical  
International



Haier Biomedical  
International



## Product Advantages



### Superior Temperature Uniformity

The temperature uniformity is  $\pm 3^{\circ}\text{C}$  with the setting temperature of  $-30^{\circ}\text{C}$



### Green and Environmental Friendly

Hydrocarbon refrigerant, energy saving and environmental protection, built with LBA foam insulation, complies with European RoHS requirements



### Microprocessor Controller

Microprocessor control, with temperature display precision at  $0.1^{\circ}\text{C}$ . Internal temperature is adjustable from  $-10^{\circ}\text{C}$  to  $-30^{\circ}\text{C}$  with  $1^{\circ}\text{C}$  control accuracy



### High-efficiency and Energy-saving Compressor and System

High-efficiency hydrocarbon compressor, optimized refrigeration system and low power consumption



### Multiple Alarm Functions

High/Low temperature/door ajar/power failure/sensor error/high ambient temperature alarm and also equipped with remote alarm interface, alarm mode is sound and flash



### Compact Footprint

Combine a small footprint with large storage capacity to provide flexible, reliable and convenient sample storage

## Optimal Design



### New Generation Appearance

A redesigned control panel features a revamped layout and updated icons, significantly improving user experience



### Low Noise Design

Noise-reduction features and structure noise reduction design



### Safe and Secure

Upgraded handle design for one-handed use, with a lockable feature for added security. Can also be paired with a padlock for multi-user convenience



### Porthole

Two portholes as standard, allows for independent testing of cabinet temperature



### Adjustable Storage Shelf Height

Built-in evaporator design, and the shelf height can be adjusted according to the size of the stored articles, so as to ensure the efficient use of large space in the cabinet



### Configuration

Standard: 5 coated steel wired shelves, remote alarm interface, USB (115V60Hz)

Optional: stainless steel shelves, blood basket, drawer, temperature recorder, RS485

**Product Parts**



**Product Dimensions**

Ramark  
W0=Total Width  
W1=Freezer Width  
W2=Inner Cabinet Width  
W3=Inner Cabinet Shelf Width  
D0=Total Depth  
D1=Inner Cabinet Width  
D2=Inner Cabinet Shelf Width  
H0=Total Height  
H1=Inner Cabinet Height

CODE MODEL	W0	W1	W2	W3	D0	D1	D2	D3	H0	H1
DW-30L698	874mm	780mm	640mm	604mm	930mm	755mm	722mm	691mm	1984mm	1458mm

 **Specifications**

Model		DW-30L698		
Technical Data	Cabinet Type	Upright		
	Climate Class	SNN		
	Cooling Type	Direct Cooling		
	Defrost Mode	Manual		
	Refrigerant	HC		
	Sound Level (dB(A))	41	43	
Performance	Cooling Performance (°C)	-30		
	Temperature Range(°C)	-10~-30		
Control	Controller	Microprocessor		
	Display	LED		
Electrical Data	Power Supply (V/HZ)	220-240/50	115/60	
	Electrical Current (A)	3.1	5	
Dimensions	Capacity (L/Cu.Ft)	698/24.7		
	Net/Gross Weight (approx)	(kg)	180/225	
		(lbs)	396.8/496	
	Interior Dimension (W*D*H)	(mm)	640*755*1458	
		(in)	25.2*29.7*57.4	
	Exterior Dimension (W*D*H)	(mm)	874*930*1984	
		(in)	34.4*36.6*78.1	
	Packing Dimension (W*D*H)	(mm)	950*1060*2110	
(in)		37.4*41.7*83.0		
	Container Load (20'/40'/40'H)	10/22/22		
Alarms	High/Low Temperature	Y		
	Remote Alarm	Y		
	Power Failure	Y		
	Sensor Error	Y		
	Low Battery	Y		
	High Ambient Temperature	Y		
	Dirty Condenser	N		
	Door Ajar	Y		
Accessories	Caster	Y		
	Leveling Feet	Y		
	Porthole	Y		
	Drawers/Inner Doors	Optional/3		
	Adjustable Shelves	Y (5)		
	USB Interface	Optional	Y	
	RS485	Optional		
	NFC	N		
Temperature Recorder	Optional			
Other	Certificate	CE	UL, ENERGY STAR	

\*Haier Biomedical reserves the right to change products and specifications without prior notice.