Your Haier Biomedical Partner



No.280 Feng Yuan Road, High-tech Zone, Qingdao, 266111, P.R. China E-mail: inquiry@haierbiomedical.com Website: www.haiermedical.com














r Biomedical Ha ternational

lote: If a slight difference occurs between pictures and actual products, please refer to actual products. Our company

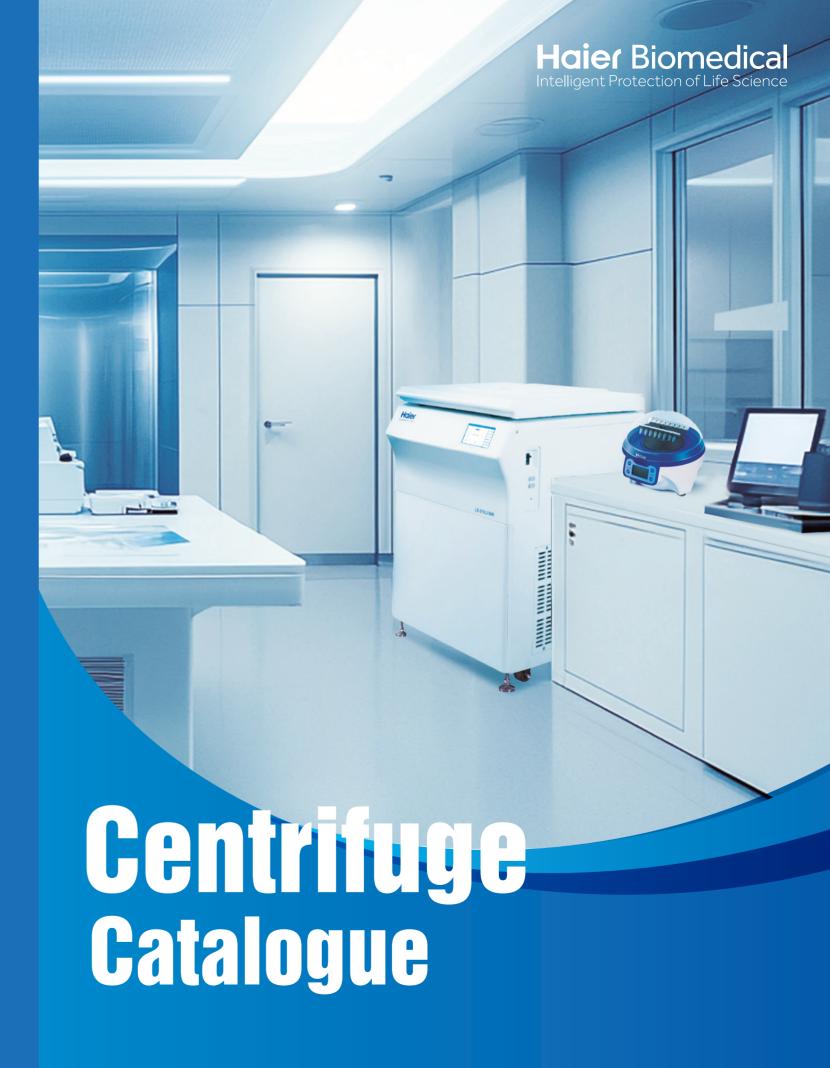



Floor-Standing Low-Speed Refrigerated Centrifuge _X-75L2400R ······	01
Floor-Standing Low-Speed Refrigerated Centrifuge  _X-60L1000R	11
Floor-Standing Highspeed Refrigerated Centrifuge _X-210L750R ······	19
Floor-Standing Large-Capacity Refrigerated Centrifuge  _X-100L1000R	33
Benchtop High-Speed Refrigerated Microcentrifuge _X-165T2R ·····	41
Benchtop High-Speed Microcentrifuge _X-165T2-J ······	45
Benchtop Multipurpose Refrigerated Centrifuge _X-155T500R······	49

Benchtop High-Speed Refrigerated Centrifuge  LX-185T100R	<b>57</b>
Benchtop High-Speed Centrifuge LX-185T100-J	- <b>67</b>
Benchtop Low-Speed Centrifuge LX-60T100-J	75
Benchtop Low-Speed Large-Capacity Centrifuge LX-60T500-J	83
Mini-Centrifuge LX-120T2Z / 50T2Z ······	93
Centrifuge Consumables	97

## Haier Biomedical Intelligent Protection of Life Science

Category	Model	Max capacity	Max speed(rpm)	Max RCF(xg)	Refrigerated	Page
Floor-Standing Blood Banking	LX-75L2400R	8*2000ml or 16*500ml blood bags	7,500	13,376	Y	01
Floor-Standing Blood Banking	LX-60L1000R	6*1000ml or 6*500ml blood bags	6,000	6,279	Y	11
Floor-Standing Bioprocessing	LX-210L750R	4*750ml	21,000	47,400	Y	19
Floor-Standing Bioprocessing	LX-100L1000R	6*1000ml	10,000	18,590	Y	33
Micro	LX-165T2R	12*10ml	16,500	25,872	Y	41
Micro	LX-165T2-J	12*5ml	16,500	25,872	N	45
Multipurpose	LX-155T500R	4*500ml	15,500	26,322	Y	49

Category	Model	Max capacity	Max speed(rpm)	Max RCF(xg)	Refrigerated	Page
Multipurpose	LX-185T100R	6*100ml	18,500	30,130	Y	57
Multipurpose	LX-185T100-J	6*100ml	18,500	30,130	N	67
Multipurpose	LX-60T100-J	6*100ml	6,000	4,749	N	75
Multipurpose	LX-60T500-J	4*500ml	6,000	5,353	N	83
Mini	LX-120T2Z	8*1.5/2ml	12,000	7,325	N	93
Mini	LX-50T2Z	8*1.5/2ml	5,000	1,270	N	93



## Floor-Standing Low-Speed Refrigerated Centrifuge









LX-75L2400R

## -50

## Advanced Centrifugation Solutions for Modern Blood Banking

Step into a realm of unmatched convenience with the Haier Biomedical LX-75L2400R Blood Banking Centrifuge, tailored to exceed the demanding needs of blood bank centrifugation. Redefine efficiency as the LX-75L2400R effortlessly accommodates up to  $16 \times 500$  ml blood bags or  $6 \times 2400$  ml centrifuge bottles, delivering optimal sedimentation results with each run

Elevate your laboratory's capabilities with the LX-75L2400R's extensive rotor options, its formidable motor, and user-friendly control panel, all complemented by a robust refrigeration system and an array of versatile accessories

Crafted for excellence, the LX-75L2400R is the centrifuge of choice for blood banks and transfusion centers. Its superior design also finds extensive application in the forefront of biopharmaceuticals, cell therapy, biological culture, vaccine production, and academic research, making it a pivotal tool in diverse scientific fields

03/04





Unlock the full potential of your blood processing center with the LX-75L2400R's expansive capacity and cutting-edge selection of rotors, carriers, and adapters, masterfully designed to handle a diverse array of high-capacity applications. This powerhouse is the epitome of versatility, perfect for sophisticated bioprocessing, clinical diagnostics, and comprehensive blood sample preparation

Quiet yet powerful, the LX-75L2400R redefines operational convenience, setting a new standard in the industry. At its heart lies a microprocessor-controlled frequency motor, delivering robust torque and rapid acceleration, all while ensuring an extended motor lifespan. Achieve peak performance with speeds soaring up to 7,500 rpm and forces reaching 13,376 x g, propelling your research to new heights with unmatched precision and reliability



## **Expansive Memory**

Maximize efficiency with over 1000 program memories, ready to streamline your workflow



## Rapid Cooling

Achieve 4°C in just 15 minutes from room temperature, speeding up your readiness



#### **Quiet Operation**

Work in peace with a sound design that keeps noise levels below 60 dB, even at maximum speed



#### **Precise Control**

Tailor your operations with 15 acceleration and 16 deceleration stages for unparalleled precision



#### Maintenance-Free Motor

Enjoy seamless operation with our brushless induction motor that requires no maintenance and ensures a clean lab environment



#### Stability

Trust in the stability of our low-center-of-gravity design for consistent performance and optimal centrifugal results



#### One-Touch Automation

Experience unparalleled convenience with our one-button automatic door mechanism that seamlessly opens and closes from front to back, ensuring quick and easy access while maintaining sterility standards.



#### Automated Rotor Lid Management

With our centrifuge, the rotor lid opens and closes automatically, complete with lid storage, reducing manual effort and enhancing workflow efficiency.



#### Intuitive 10-Inch Touch Screen

The expansive 10-inch touch screen interface provides operators with user-friendly navigation through settings and features, making operation straightforward and minimizing training time.



#### **Electric Push Rod**

Engineered for safety, the electric push rod offer smooth and controlled opening and closing of the centrifuge chamber, providing an extra layer of protection to users.



#### 316 Stainless Steel Chamber

The premium 316 stainless steel construction of the centrifuge chamber ensures a rust-free and contaminant-resistant environment that is exceptionally easy to clean and maintain.



#### **Versatile Rotor Capacity**

Our centrifuge accommodates a variety of research needs with three rotors designed for 250ml, 500ml, 1000ml, 2000ml, and 2400ml centrifuge bottles, as well as 500ml blood bags, maximizing your laboratory's throughput capabilities.



#### **Emergency Lid Release**

The lid can be opened by utilizing the reset hole on the side plate. Rotate a 5mm hex key clockwise several times to open the lid.



## **Interactive 10-Inch LCD Touchscreen**



## Glove-Friendly and Ergonomic

The 45° angled, glove-friendly touchscreen offers a wide viewing angle, making it comfortable and practical for users of all heights



## **Enhanced Access Controls**

Equipped with user login and password protection, ensuring individualized program data management for each user



Multiple operation setting modes, including unique blood bank mode function, convenient for blood bank users to use



## **Visibility and Monitoring**

Enlarged track bar display of set and actual run conditions during centrifugation, ensuring clear visibility from a distance



## **Program Recall and Efficiency**

Quick recall feature for over 1000+ programs with alpha-numeric naming, minimizing time between runs and facilitating quick starts for routine users



## Adjustable On-the-Fly

Users can alter speed, centrifugal force, time, and temperature during the centrifugation process for precise control







Multi-echelon program running mode, achieving continuous operation of multiple speed steps for special users



Electronic level assists in installation and leveling, making installation easier and simpler

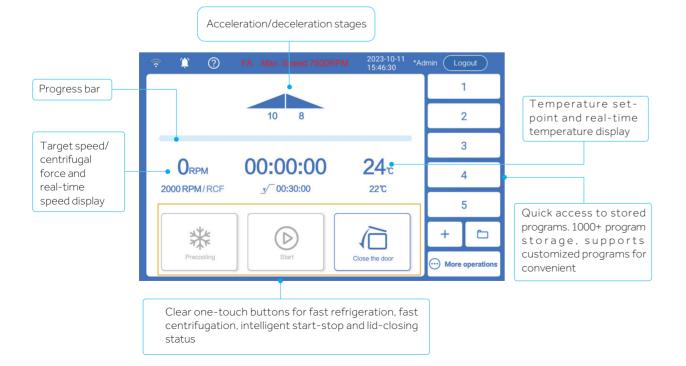

The operation record including speed, time and temperature, operation curve with alarm messages etc. can be tracked on the screen and exported by **USB** port

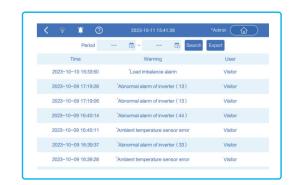







#### Immediate Imbalance Detection & Shutdown

Offers real-time loading imbalance detection with instant shutdown to protect against imbalances



#### Stable Installation

 $\Phi 120 mm$  adjustable feet ensure the centrifuge is securely and stably installed in any setting



## **Comprehensive Alarm System**

Features a multi-faceted alarm system for abnormal voltage, high temperature, overspeed, load imbalance, rotor absence, brake anomalies, and rotor life



#### **Anti-Vibration Lid Lock**

Prevents accidental opening during operation, ensuring sample integrity and operator safety



### **Robust Construction**

Fully covered steel plate protective sleeve and a 12mm steel shell provide enhanced durability and user safety




## **High-Quality Connectivity**

Equipped with an imported 5-core aviation plug for easy, fast, and safe usage




## Blood Station User: Blood centrifugal test

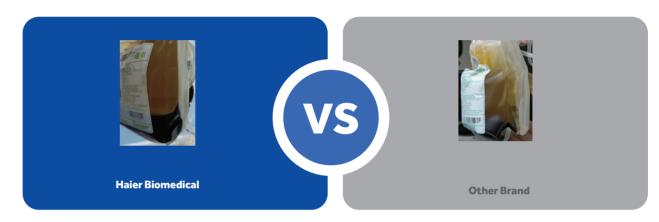

## **Precision Layering**

Clear Separation with Distinct Leukocyte Layer, Eliminating the Need for Re-Separation

## Scientific Research User: Seaweed centrifugal test

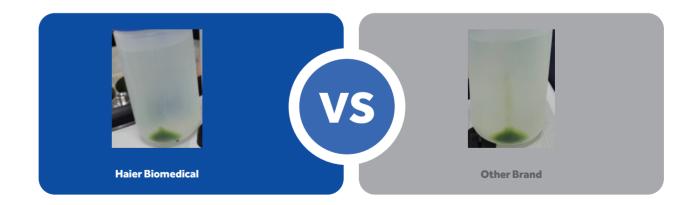

### Optimal Centrifugation

Clear Layering, High Yield, and Zero Wall Adhesion



## Technical Data

Model	LX-75L2400R
Maximum Capacity	8*2000ml* or 16*500ml blood bags
Max Speed (rpm)	7,500
Max RCF(xg)	13,376
Speed Accuracy (rpm)	±5
Time Setting Range	1min~99h59m59s (1-second increment)
Precooling Function	Yes
Temperature Setting Range (°C)	-20~40, adjustment increment: 1
Display	10 inch touch screen
Noise (dB(A))	≤60
Power Supply	3~400V 50/60Hz
Power (W)	6400
Exterior Dimensions (W*D*H)(mm)	900*1050*1020
Lid Opening Height (mm)	1780
Net Weight (kg)	620
Packing Dimensions (W*D*H)(mm)	1092*1279*1720
Gross Weight (kg)	750
Other Functions	Rotation speed/RCF switchover, inching operation, running process indicator, sound alert, adjustable acceleration/deceleration, user logging, user lock-out, automatic door opening, on-screen display for imbalance, over temperature, stainless steel chamber
Certification	CE

<sup>\*</sup>Buckets for 2000ml centrifuge bottles can be customized depending on requirements.











BE No.	BE12DC000	BE12JP002	BE12JQ002	BE12K1000
Туре	Fixed-angle rotor (Sealed)	Adapter (Sets of 6) for BE12DC000	Adapter (Sets of 6) for BE12DC000	Centrifuge Bottle(Sets of 6) for BE12DC000
Max Capacity	6*1000ml	500ml	250ml	1000ml
Max RCF (xg)	13,376	-	-	-
Max Speed (rpm)	7,500	-	-	-
Max Bottle Capacity (ml)	1000	500	250	-
Bottles Per Rotor	6	6	6	-
Bottle Type	Flat bottom	Flat bottom	Flat bottom	-

## Swing-Out Rotor



Swing-out Rotor 12 blood bags

Bucket (Sets of 6)



BE12H4000

2\*500ml blood bags

or 2400ml

Flat bottom

BE No.	BE12F5001		BE No.
Max. RCF (xg)	5,331	4	Max Capacity
Max. speed (rpm)	4,000		Bucket Type





#### Adapter/Bottle(Sets of 6) for BE12H4000

BE No.	BE12JT001	BE12K0000	
Туре	Double blood bag adapter	Centrifuge bottle	
Max Capacity	2*500ml blood bags	2400ml	
Blood bags/Bottles per rotor	12 blood bags	6 bottles	













BE No.	BE12H9000
Max Capacity	2*500ml blood bags
Bucket Type	Flat bottom

#### Adapter(Sets of 8) for BE12H9000

- 1. 1. 1. 1. 1. 1. 1. 1. 1. 1. 1. 1. 1.			
BE No.	BE12JF002		
Туре	Double blood bag adapter		
Max Capacity	2*500ml blood bags		
Blood bags per rotor	16		







## Floor-Standing Low-Speed Refrigerated Centrifuge

## **Scope of Application**

Haier Biomedical LX-60L1000R Blood Banking Centrifuge seamlessly integrates a remarkably compact footprint with top-tier performance. It can handle up to 6×500ml blood bags at a time, providing excellent accuracy in speed and temperature control. This ensures optimal sample integrity and consistency throughout each centrifugation cycle.



## **Excellent Performance**



#### **Precision Control**

Ensure excellent speed and temperature control accuracy, and maintain the integrity and consistency of samples in each centrifugation cycle, thus achieving high-quality blood treatment.



### **Expansive Memory**

Maximize efficiency with over 1000 program memories, enabling the streamlining of the workflow.



### **Rapid Cooling**

The unique refrigeration technology enables rapid cooling. It achieves 4°C in just 10 minutes from room temperature under pre-cooling mode, offering impressive time savings.



### **Automated Door Management**

With one-touch automatic door opening as a standard feature, optional automatic door opening/closing functionality is also available



## **Quick-Lock Rotor Exchange**

Patented tool-free rotor exchange technology enables the installation and removal of rotors on the centrifuge automatically and securely without the need for tools.



## **Quiet Operation**

The sound design is optimized for low noise levels ( $\leq$  54dB(A)).



## Ergonomic Design



A space-efficient rotor lid storage design



Optional electric push rod offer smooth and controlled opening and closing of the centrifuge door



7-inch touch screen interface for convenient operation.



Quick rotor exchange design for installation and removal with no tools required.



316 stainless steel centrifuge chamber is rust-free and easy to clean.



USB port for exporting operation records, including data on speed, time, temperature, operation curves, and alarm messages.



Emergency lid release: Pry out the side - plate's white plug with a tool and pull the release cord 7+ times to open the lid.



## Large, bright, interactive 7-inch LCD interface offers:

The 45 ° angled, glove-compatible touchscreen provides a wide viewing angle for easy visibility and comfortable use by technicians of all heights.

Advanced user settings and password protection ensure safe operation. Multiple user profiles can be created to customize access and settings, improving workflow and sample safety.

During operation, the display of settings and actual operating conditions can be easily monitored through the enlarged tracking bar, so as to obtain a clear vision and ensure that the operator can see them from a distance.

Store more than 1,000 programs, and can directly access 5 or 10 commonly used programs, depending on the model. Programs can be named alphanumerically, ensuring easy identification and selection, minimizing the risk of errors during operation.

Users can change the speed, centrifugal force, time and temperature during the centrifugation.





Multiple operation setting modes, including unique blood bank mode function, convenient for blood bank users to use.




Electronic level assists in installation and leveling, making installation easier and simpler.



Multi-echelon program running mode, achieving continuous operation of multiple speed steps for special users



The operation record including speed, time and temperature, operation curve, alarm messages etc. can be tracked on the screen and exported by USB port.

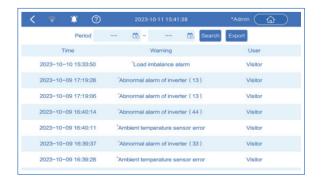







Automatic radio - frequency rotor identification (RFID) enables immediate identification of a rotor once it is securely placed in the centrifuge chamber, avoids over-speed, reduces error messages, and enhances the safety.



Real time loading imbalance detection and intelligent lid locking monitoring, immediate shutdown protection in case of imbalance or lid opening.



Use of brushless induction motors to provide maintenance-free operation, eliminating carbon particle emissions.



Multiple alarms including abnormal voltage alarm, high temperature alarm, overspeed alarm, load imbalance alarm, rotor absent alarm, brake abnormal alarm, rotor life reminder alarm, and other alarm functions can effectively ensure the safety of the centrifuge.



A triple - layer construction paired with a fully covered steel plate protective sleeve provides an effective safeguard.



Model	LX-60L1000R
Max Capacity	6*1000ml or 6*500ml blood bags
Speed Rang (rpm)	500-6,000
Max RCF (xg)	6,279
Speed Accuracy (rpm)	±5
Time Setting Range	$1$ min $\sim$ 99h59m59s (1-second increment)
Precooling Function	Yes
Temperature Setting Range	-20°C~40°C, adjustable to 1°C
Display	7 inch touch screen
Noise (dB (A))	≤54
Power Supply (V/Hz)	220~240/50/60
Power (W)	3500
Exterior Dimensions (W*D*H) (mm)	740*920*990
Lid Opening Height (mm)	1720
Net Weight (kg)	338
Packing Dimensions (W*D*H) (mm)	940*1150*1270
Gross Weight (kg)	423
Certification	CE





BE No.	BE12DL001
Туре	Sealed
Max Capacity	6*500ml
Max RCF (xg)	6,279
Max Speed (rpm)	6,000















BE12DL001										
BE No.	No adapter needed	BE12JT002	BE12JU002	BE12JV002	BE12JW002	BE12JX002	BE12JY002	BE12JZ002	BE12J0003	BE12J1003
Max Capacity	500ml	250ml	100ml	2*50ml	2*50ml	4*30ml	7*15ml	8*10ml	10*5ml	14*1.5/2ml
Tubes Per Adapter/Rotor	r 6	1/6	1/6	2/12	2/12	4/24	7/42	8/48	10/60	14/84
Tube					Опининин					
Max Tube Dimensions (Ø x L mm)	69.5*170	62*123	38*120	29*105	29*118	25*102	17*121	16*86	13*54	11*38
Tube Type	Flat bottom	Flat bottom	Round bottom	Round bottom	Conical bottom	Round bottom	Conical bottom	Round bottom	Round bottom	Conical/Round bottom











BE No.	BE12FT001	
Max. RCF (xg)	5,100	4
Max. speed (rpm)	4.200	





BE No.	BE12G9000
Max Capacity (ml)	1000
<b>Bucket Type</b>	Flat bottom

## Adapter(Sets of 6) for BE12G9000













BE No.	No adapter needed	BE12J8004	BE12J9004	BE12J7004	BE12J6004	BE12J5004	BE12JA004
Max Capacity	1000ml	500ml	500ml	300ml	250ml	200ml	4*100ml
Bottles Per Adapter /Rotor	- 6	1/6	1/6	1/6	1/6	1/6	4/24
Tube							
Max Tube Dimensions (Ø x L mm)	102*179	96*150	80*136	60*138	60*165	/	38*120
Tube Type	Flat Bottom	Conical Bottom	Flat Bottom	Flat Bottom	Conical Bottom	Blood Bag	Round Bottom

















BE No.	BE12JB004	BE12JC004	BE12JD004	BE12JE004	BE12JH004	BE12JF004	BE12JG004
Max Capacity	5*50ml	5*50ml	14*15ml	14*15ml	37*7ml	24*5/7ml	28*1.5/2ml
Tubes Per Adapter /Rotor	5/30	5/30	14/84	14/84	37/222	24/144	28/168
Tube	Management (						
Max Tube Dimensions (Ø x L mm)	29*118	29*105	17*121	17*121	12*100	13*75/106	11*38
Tube Type	Conical Bottom	Round Bottom	Conical Bottom	Round Bottom	Round Bottom	Small/Large Blood Collection Tube	Conical / Round Bottom











BE No.	BE12HD000
Max Capacity (ml)	1000
Bucket Type	Flat bottom

## Adapter(Sets of 6) for BE12HD000

BE No.	BE12J3004	
Туре	Blood bag adapter	
Max Capacity	500ml blood bag	
Blood Bags Per Adapter /Rotor	1/6	









LX-210L750R

# Floor-Standing Highspeed Refrigerated Centrifuge

## **Scope of Application**

Haier Biomedical LX-210L750R is a floor-standing highspeed refrigerated centrifuge features with 3L capacity and performance up to 47,400 x g, which is ideal for all-purpose use in the industrial and research sectors.





Offering a broad portfolio of fixed-angle and swing-out rotors, which can accommodate tubes and bottles from 1.5 ml to 750 ml, as well as deep-well plates



Automatic radio-frequency rotor identification (RFID) that enables immediate identification of a rotor once it is securely placed in the centrifuge chamber, avoid over-speed, reduce error messages, and enhance the safety



Low center of gravity design, amplitude less than 0.02mm, ensuring stable operation and good centrifugal effect



Use of brushless induction motors to provide maintenance-free operation, eliminating carbon particle emissions



High-strength, aero-aluminum rotor, which is acid and alkali corrosion-resistant, greatly increases the service life, and allows for 50,000 uses



Optimized multi-layer shock absorption technology makes the sound level ≤ 56(dB(A))



## Product Parts



Dual self-lock design to ensure safe operation



Dual gas springs with safety protection



4 universal casters and 4 levelling feet allows for easy movement and locks safely into place



316 stainless steel centrifuge chamber is rust-free and easy to clean

USB port for exporting operation record including speed, time and temperature, operation curve with alarm messages etc



• Patented tool-free rotor exchange technology enables the installation and ensure safe use



technology enables the installation and removal of rotors on the centrifuge automatically and securely without the need for an tools





## Intuitive Interface



The 7-inch 45 ° glove-friendly touchscreen, the wide viewing angle more comfortable and practical for users of all heights



Easily monitor the display of set and actual run conditions with enlarged track bar during the run for clear visibility, which ensures the operator can see them from a distance



The speed, centrifugal force, time and temperature can be changed by user during the centrifugation process



Further enhanced with access controls such as user login with password protection, the independent program data of each user is managed separately



Quick recall of up to 1000+ programs with alpha-numeric naming to minimize time between runs and for quick run start by routine users



11 acceleration + 12 deceleration stages



fast centrifugation, intelligent start-stop and lid-closing status





Real time loading imbalance detection and intelligent lid locking monitoring, immediate shutdown protection in case of imbalance or lid opening



Multiple alarm including abnormal voltage alarm, high temperature alarm, overspeed alarm, load imbalance alarm, rotor absent alarm, brake abnormal alarm, rotor life reminder alarm, and other alarm functions



Adopting  $\Phi$ 80mm adjustable feet for more stable installation



Fully covered steel plate protective sleeve and 12mm steel shell



Model	LX-210L750R
Max Capacity	4*750ml
Speed Range (rpm)	500-21,000
Max RCF (xg)	47,400
Speed Accuracy (rpm)	±10
Time Setting Range	1min~99h59m59s(1-second increment)
Precooling Function	Yes
Temperature Setting Range (°C)	-20~40, adjustable to 1
Display	7 inch touch screen
Noise (dB(A))	≤56
Power Supply (V/Hz)	220~240/50/60
Power (W)	3400
Exterior Dimensions (W*D*H)(mm)	740*920*990
Lid Opening Height (mm)	1720
Net Weight (kg)	375
Packing Dimensions (W*D*H)(mm)	940*1150*1270
Gross Weight (kg)	470
Certification	CE







BE No.	BE12DM001
Туре	Carbon fiber, Sealed
Max. Capacity (ml)	4*1000
Max RCF (xg)	16,902
Max Speed (rpm)	9,400









## Adapter(Sets of 6) for BF12DM001

Adapter(Sets of 6) for BE	12DM001				
BE No.	No adapter needed	BE12JP002	BE12JQ002	BE12JR002	BE12JS002
Max Capacity	1000ml	500ml	250ml	6*50ml	5*50ml
Tubes Per Adapter/Rotor	4	1/4	1/4	6/24	5/20
Tube	*		Ō		
Max Tube Dimensions (Ø x L mm)	97.5*195	69.5*170	62*123	29*105	29*118
Tube Type	Flat bottom	Flat bottom	Flat bottom	Round bottom	Conical bottom

\*Haier Biomedical offers 1000ml centrifuge bottles(sets of 6) - BE12K1000







BE No.	BE12D0001
Туре	Sealed
Max Capacity	6*500ml
Max RCF (xg)	17,778
Max Speed (rpm)	10,000

















BE No.	No adapter	PE12 IT002	PE12 II I002	PE12 IV002	BE12JW002	RE12 IV002	RE12 IV002	BE12JZ002	BE12J0003	BE12J1003
Max Capacity	needed 500ml	250ml	100ml	2*50ml	2*50ml	4*30ml	7*15ml	8*10ml	10*5ml	14*1.5/2ml
Tubes Per Adapter/Roto	r <sup>6</sup>	1/6	1/6	2/12	2/12	4/24	7/42	8/48	10/60	14/84
Tube							-Cilmononono			
Max Tube Dimensions (Ø x L mm)	69.5*170	62*123	38*120	29*105	29*120	25*102	17*121	16*82	13*54	11*38
Tube Type	Flat bottom	Flat bottom	Round bottom	Round bottom	Conical bottom	Round bottom	Conical bottom	Round bottom	Round bottom	Conical /Round bottom







BE No.	BE12D1001
Туре	Sealed
Max Capacity	6*250ml
Max RCF (xg)	25,115
Max Speed (rpm)	12,000



BE No.	BE12D5001
Туре	Sealed
Max Capacity	4*250ml
Max RCF (xg)	29,440
Max Speed (rpm)	14,000











BE No.	No adapter needed	BE12JF003	BE12JG003	BE12JH003	BE12JJ003
Max Capacity	250ml	100ml	50ml	50ml	5*15ml
Tubes per adapter/ rotor (BE12D1001)	6	1/6	1/6	1/6	5/30
Tubes per adapter/ rotor (BE12D5001)	4	1/4	1/4	1/4	5/20
Tube					
Max Tube Dimensions (Ø x L mm)	62*123	38*120	29*105	29*120	17*121
Tube Type	Flat bottom	Round bottom	Round bottom	Conical bottom	Conical botton





BE No.	BE12D4001
Туре	Sealed
Max Capacity	10*100ml
Max RCF (xg)	21,280
Max Speed (rpm)	12,000

## Adapter(Sets of 10) for BE12D4001











BE No.	No adapter needed	BE12JS003	BE12JR003	BE12JQ003	BE12JN003	BE12JP003	BE12JT003
Max Capacity	100ml	50ml	50ml	25ml	15ml	15ml	2*10ml
Tubes Per Adapter/Rotor	. 10	1/10	1/10	1/10	1/10	1/10	2/20
Tube							
Max Tube Dimensions (Ø x L mm)	38*120	29*105	29*120	22*110	17*121	17*121	16*82
Tube Type	Round bottom	Round bottom	Conical bottom	Round bottom	Round bottom	Conical bottom	Round bottom







BE No.	BE12D7001
Туре	Sealed
Max Capacity	8*50ml
Max RCF (xg)	38,460
Max Speed (rpm)	18,000



#### Adapter(Sets of 8) for BE12D7001

BE No.	No adapter needed	BE12JK003
Max Capacity (ml)	50	15
Tubes Per Adapter/Rotor	1/8	1/8
Tube		
Max Tube Dimensions(Ø x L mm)	29*105	17*121
Tube Type	Round bottom	Round bottom



BE No.	BE12D6001
Туре	Sealed, round bottom
Max Capacity	6*50ml
Max RCF (xg)	44,800
Max Speed (rpm)	20,000

## Adapter(Sets of 6) for BE12D6001





BE No.	No adapter needed	BE12JM003	BE12JL003
Max Capacity (ml)	50	30	15
Tubes Per Adapter/Rotor	1/6	1/6	1/6
Tube			
Max Tube Dimensions (Ø x L mm)	29*105	25*102	17*121
Tube Type	Round bottom	Round bottom	Round bottom





BE No.	BE12D8001
Туре	Sealed, conical bottom
Max Capacity	6*50ml
Max RCF (xg)	44.800
Max Speed (rpm)	20,000

## Adapter(Sets of 6) for BE12D8001



BE No.	No adapter needed	BE12JC002
Max Capacity (ml)	50	15
Tubes Per Adapter/Rotor	1/6	1/6
Tube		
Max Tube Dimensions (Ø x L mm)	29*120	17*121
Tube Type	Conical bottom	Conical bottom







BE No.	BE12D3001
Туре	Sealed
Max Capacity	12*15ml
Max RCF (xg)	28,618
Max Speed (rpm)	16,000

Tube for BE12D3001	The state of the s
Max Tube Capacity (ml)	15
Tubes per Rotor	12
Max Tube Dimensions (Ø x L mm)	17*121
Tube Type	Conical bottom



BE No.	BE12D2001
Туре	Sealed
Max Capacity	16*10ml
Max RCF (xg)	47,400
Max Speed (rpm)	21,000

Tube for BE12D2001	
Max Tube Capacity (ml)	10
Tubes per Rotor	16
Max Tube Dimensions (Ø x L mm)	16*82
Tube Type	Round bottom











BE No.	BE12FS001	
Max. RCF (xg)	3,488	7
Max. speed (rpm)	4,000	

BE No.	BE12G5000
Max. RCF (xg)	3,488
Max. speed (rpm)	4,000
Max Capacity (ml)	750
Bucket Type	Flat bottom















BE No.	No adapter needed	BE12J3003	BE12J4003	BE12J5003	BE12J6003	BE12J7003	BE12J8003	BE12J9003	BE12JA003
Max Capacity	750ml	500ml	250ml	4*100ml	7*50ml	16*15ml	23*10ml	31*5/7ml	56*1.5/2ml
Tubes Per Adapter/Rotor	4	1/4	1/4	4/16	7/28	16/64	23/92	31/124	56/224
Tube					Финиппини	Tipopopopopopo			
Max Tube Dimensions (Ø x L mm)	99*140	96*150	60*165	38*120	29*118	17*121	16*106	13*75/106	11*38
Tube Type	Flat bottom	Conical bottom	Conical bottom	Round bottom	Conical bottom	Conical bottom	Large blood collection tube ( <b>Φ</b> 16mm)	Short/long blood collection tube	Conical/ Round bottom















BE No.	BE12FS001
Max. RCF (xg)	3,488
Max. speed (rpm)	4,000

BE No.	BE12HB000
Max. RCF (xg)	3,488
Max. speed (rpm)	4,000
Max Capacity	-
Bucket Type	Flat bottom









Adapter(Sets of 4) for BE12HB000				
BE No.	BE12JB003	BE12JC003	BE12JD003	BE12JE003
Max Capacity	9*50ml	20*15ml	30*10ml	39*5/7ml
Tubes Per Adapter/Rotor	9/36	20/80	30/120	39/156
Tube	(http://doi.org/	Distribution (2001).		
Maximum Tube Dimensions (Ø x L mm)	29*118	17*121	16*106	13*75/106
Tube Type	Conical bottom	Conical bottom	Large blood collection tube ( <b>Φ</b> 16mm)	Short/long blood collection tube












CIRCLE CONTRACTOR	BE No.	BE12FS001
	Max. RCF (xg)	3,488
0	Max. speed (rpm)	4,000

BE No.	BE12HC000
Max. RCF (xg)	3,148
Max. speed (rpm)	4,000
Max Capacity	-
Bucket Type	Rectangular bottom Max height 60mm







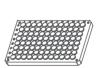





Plates per bucket/rotor	4/16	3/12	3/12	1/4	1/4
W*D*H (mm)	128*86*15	128*86*17.5	128*86*20	128*86*44.5	128*86*60
Plate Type	MTP (without lid)	MTP (with lid)	PCR plate 96-fold	DWP	Filter plate





# Floor-Standing Large-Capacity Refrigerated Centrifuge

## **Scope of Application**

Haier Biomedical LX-100L1000R is a conventional centrifuge designed for biopharmaceuticals with large capacity 6L, which is widely used in fields such as biopharmaceuticals and academic research. It plays an important role in the processes of protein collection, separation, and purification for samples like macromolecular proteins and bacteria.





Satisfy the changing application requirements that extend from academic research to production facilities. Offer up to  $18,590 \times g$  of centrifugal force and a maximum capacity of 6 liters



Automatic radio-frequency rotor identification (RFID) that enables immediate identification of a rotor once it is securely placed in the centrifuge chamber, avoid over-speed, reduce error messages, and enhance the safety



Optimized multi-layer shock absorption technology makes the sound level ≤ 63 (dB (A))



Low center of gravity design, amplitude less than 0.02mm, ensuring stable operation and good centrifugal effect



## Product Parts



Dual self-lock design to ensure safe operation



Dual gas springs with safety protection



4 universal casters and 4 levelling feet allows for easy movement and locks safely into place



316 stainless steel centrifuge chamber is rust-free and easy to clean



USB port for exporting operation record including speed, time and temperature, operation curve with alarm messages etc



Triple-layer construction to ensure safe use



Emergency lid release, pry out the white plastic plug on the side plate with a tool and pull the attached release cord to open the lid



Offering 6\*1000ml, 6\*500ml fixed-angel rotor with carbon fiber/aluminum alloy for option, and a variety of adapters available, which can accommodate tubes and bottles from 1.5 ml to 1 L



- The high-efficiency carbon fiber rotor achieves the maximum centrifugal force, faster acceleration and deceleration times, and greatly improves centrifugal efficiency
- Corrosion resistance and high strength can ensure the security and durability of the rotor
- The service life can reach up to 15 years



- High-strength, aero-aluminum rotor, which is acid and alkali corrosion-resistant, greatly increases the service life
- A comprehensive adapter portfolio, which is compatible with 500ml, 250ml, 100ml, 50ml, 15ml, 10ml, 5ml tubes and bottles
- The safe usage times can reach up to 50,000





The 7-inch 45° glove-friendly touchscreen, the wide viewing angle more comfortable and practical for users of all heights



Easily monitor the display of set and actual run conditions with enlarged track bar during the run for clear visibility, which ensures the operator can see them from a distance



The speed, centrifugal force, time and temperature can be changed by user during the centrifugation process



Further enhanced with access controls such as user login with password protection, the independent program data of each user is managed separately



Quick recall of up to 1000+ programs with alpha-numeric naming to minimize time between runs and for quick run start by routine users  $\frac{1}{2}$ 



11 acceleration + 12 deceleration stages









Real time loading imbalance detection and intelligent lid locking monitoring, immediate shutdown protection in case of imbalance or lid opening



Multiple alarms including abnormal voltage alarm, high temperature alarm, overspeed alarm, load imbalance alarm, rotor absent alarm, brake abnormal alarm, rotor life reminder alarm, and other alarm functions



Adopting  $\Phi 80$ mm adjustable foot for more stable installation



Fully covered steel plate protective sleeve and 12mm steel shell



Use of brushless induction motors to provide maintenance-free operation, eliminating carbon particle emissions



## **Technical Data**

Model	LX-100L1000R			
Max Capacity	6*1000ml			
Speed Range (rpm)	500-10,000			
Max RCF (xg)	18,590			
Speed Accuracy (rpm)	±10			
Time Setting Range	1min-99h59m59s (1-second increment)			
Precooling Function	Yes			
Temperature Setting Range (°C)	-20~40, adjustable to 1			
Display	7 inch touch screen			
Noise (dB (A))	≤63			
Power Supply (V/Hz)	400V 3N~50/60			
Power (W)	8800			
Exterior Dimensions (W*D*H) (mm)	740*920*990			
Lid Opening Height (mm)	1720			
Net Weight (kg)	380			
Packing Dimensions (W*D*H) (mm)	940*1150*1270			
Gross Weight (kg)	475			
Certification	CE			

Laboratory centrifuges contain fluorinated greenhouse gases in the sealing equipment.





BE No.	BE12DW000
Туре	Sealed
Max Capacity	6*1000ml
Max RCF (xg)	14,472
Max Speed (rpm)	8,000



Adapter (Sets of 6) for BE12DW000, BE12DY000

BE No.	BE12DY000
Туре	Carbon fiber, Sealed
Max Capacity	6*1000ml
Max RCF (xg)	16,660
Max Speed (rpm)	8,600









BE No.	No adapter needed	BE12JP002	BE12JQ002	BE12JR002	BE12JS002
Max Capacity	1000ml	500ml	250ml	6*50ml	5*50ml
Tubes per adapter/rotor	6	1/6	1/6	6/36	5/30
Tube	*				
Max Tube Dimensions (Ø x L mm)	97.5*196	69.5*170	62*123	29*105	29*120
Tube Type	Flat bottom	Flat bottom	Flat bottom	Round bottom	Conical bottom

\*Haier Biomedical offers 1000ml centrifuge bottles(sets of 6) - BE12K1000







BE No.	BE12DX000
Туре	Sealed
Max Capacity (ml)	6*500
Max RCF (xg)	17,700
Max Speed (rpm)	10,000



BE No.	BE12DZ000		
Туре	Carbon fiber, Sealed		
Max Capacity	6*500ml		
Max RCF (xg)	18,590		
Max Speed (rpm)	10,000		











BE No.	No adapter needed	BE12JT002	BE12JU002	BE12JV002	BE12JW002
Max Capacity	500ml	250ml	100ml	2*50ml	2*50ml
Tubes per adapter/rotor	6	1/6	1/6	2/12	2/12
Tube					
Max Tube Dimensions (Ø x L mm)	69.5*170	62*123	38*120	29*105	29*120
Tube Type	Flat bottom	Flat bottom	Round bottom	Round bottom	Conical bottor

Adapter (Sets of 6) for BE12DX000, BE12DZ000











BE No.	BE12JX002	BE12JY002	BE12JZ002	BE12J0003	BE12J1003
Max Capacity	4*30ml	7*15ml	8*10ml	10*5ml	14*1.5/2ml
Tubes per adapter/rotor	4/24	7/42	8/48	10/60	14/84
Tube					
Max Tube Dimensions (Ø x L mm)	25*102	17*121	16*82	13*54	11*38
Tube Type	Round bottom	Conical bottom	Round bottom	Round bottom	Conical/Round bottom



## LX-165T2R

## **Benchtop High-Speed Refrigerated Microcentrifuge**

## Scope of Application

It is widely used in fields such as life sciences, biology, medical research, pharmaceuticals, geology, agriculture and forestry, petrochemicals, light industry, and metal particles, for rapid separation and extraction of trace

It also serves as IVD lab equipment for separating human blood, used in diagnostic tests like hematological, immunological screening, and cardiovascular assessment.



## **Product Features**



Short pre-cooling time of only 5minutes from room temperature to 4°C



Fast sample processing at maximum speeds of up to 16,500 rpm (25,872 x g)

Use of brushless induction motors to



Precise temperature and centrifugal speed control ensures optimal separation results

Optimized sound design for low noise



provide maintenance-free operation, eliminating carbon particle emissions



Real time loading imbalance detection, immediate shutdown protection in case of imbalance



Automatic radio-frequency rotor identification (RFID) and rotor life monitoring for maximum operational safety



## Intuitive Interface

levels (≤52dB(A))



The 7-inch 45° glove-friendly touchscreen, the wide viewing angle more comfortable and practical for users of all heights



Further enhanced with access controls such as user login with password protection, the independent program data of each user is managed separately



Quick recall of up to 1000+ programs with alpha-numeric naming to minimize time between runs and for quick run start by routine users



The speed, centrifugal force, time and temperature can be changed by user during the centrifugation process



Haier

LX-165T2R

05

Haier Biomedical Intelligent Protection of Life Science



2 acceleration + 3 deceleration stages

Micro



## **Product Parts**



Intelligent lid locking monitoring and automatic self-locking to ensure safe operation



Viewing port in the lid for observing the operating condition of the rotor



Drainage tray to prevent internal water accumulation and improve corrosion resistance



Triple-layer construction to ensure safe use



316 stainless steel centrifuge chamber is rust-free and easy to clean



Emergency lid release, pry out the white plastic plug on the side plate with a tool and pull the attached release cord at least 7 times to open the lid



USB port for exporting operation record including speed, time and temperature, operation curve with alarm messages



Model	Maximum Capacity	Speed Range(rpm)	Max RCF(xg)	Speed Accuracy (rpm)	Time Setting Range	Precooling Function	Temperature Setting Range(°C)	Display
LX-165T2R	12*10 ml	500-16,500	25,872	±10	1s ~ 99h59m59s (1-second increment)	Yes	-20 ~40, adjustable to 1	7 inch touch screen

Noise Level dB(A)	Power Supply (V/Hz)	Power(W)	Exterior Dimensions (W*D*H)(mm)	Net Weight (kg)	Packing Dimensions (W*D*H)(mm)	Gross Weight(kg)	Certifications
≤52	220~240/50/60	480	340*625*248	34	425*795*430	48	CE













BE No.	BE12C0Z00	BE11D0Z00	BE12D0Z00	BE12E0Z00	BE12DR000
Туре	Plastic cover	Sealed	Sealed	Sealed	Sealed
Max Capacity	24*1.5/2ml	24*1.5/2ml	12*5ml	4*8*0.2ml	12*10ml
Max RCF (xg)	25,872	21,382	21,407	18,891	2,753
Max Speed (rpm)	16,500	15,000	15,000	15,000	5,000
Tubes per Rotor	24	24	12	4	12
Tube					
Max Tube Capacity (ml)	1.5/2	1.5/2	5	0.2 ml 8-strip PCR	10
Max Tube Dimensions (Ø x L mm)	11*38	11*38	13.5*75	6*18	15*106
Tube Type	Round /Conical bottom	Round / Conical bottom	Round bottom	PCR	Large blood collection tube ( <b>Φ</b> 15mm)

#### Adapter(Sets of 50) for BE12DR000

BE No.	BE12JN001	
Max Capacity (ml)	5	
Tube		
Max Tube Dimensions (Ø x L mm)	13*75	
Tube Type	Short blood collection tube	

\*Haier Biomedical also can provide 0.2ml adapter (sets of 24) (BE12J2003),0.5ml adapter (sets of 24) (BE12JN002) for 1.5/2ml tube rotor

## **Benchtop High-Speed** Microcentrifuge

## **Scope of Application**

It is widely used in fields such as life sciences, biology, medical research, pharmaceuticals, geology, agriculture and forestry, petrochemicals, light industry, and metal particles, for rapid separation and extraction of trace samples. It also serves as IVD lab equipment for separating human blood, used in diagnostic tests like hematological, immunological screening, and cardiovascular assessment.



LX-165T2-J





Fast sample processing at maximum speeds of up to 16,500rpm (25,872xg)



Compact footprint with a wide array of microcentrifuge tubes from 0.2 to 5 ml



Use of brushless induction motors to provide maintenance-free operation, eliminating carbon particle emissions



Automatic rotor identification with over speed protection



Optimized sound design for low noise levels (≤58dB(A))



Real time loading imbalance detection, immediate shutdown protection in case of imbalance



Intuitive Interface



LCD keypad and control knob for easy operation, the wide viewing angle more comfortable and practical for users of all heights



Quick recall of up to 99 programs for quick run start by routine users



2 acceleration + 3 deceleration stages



The speed, centrifugal force, and time can be changed during the centrifugation process



LX-165T2-J

LX-165T2-J

**Haier** Biomedical Intelligent Protection of Life Science

Micro





## Product Parts



Patented anti-seismic door lock. In case of abnormal vibration, the cover can be locked reliably



Triple-layer construction to ensure safe use





Unique air duct structure design, the heat inside the centrifugal chamber can dissipate quickly



316 stainless steel centrifuge chamber is rust-free and easy to clean



Emergency lid release, pry out the white plastic plug on the side plate with a tool and pull the attached release cord to open the lid



Model	Maximum Capacity	Speed Range (rpm)	Max RCF(xg)	Speed Accuracy (rpm)	Time Setting Range	Display
LX-165T2-J	12*5 ml	500~16,500	25,872	±10	1s~99h/59min (1s/1min increment)	LCD

Noise Level dB(A)	Power Supply (V/Hz)	Power(W)	Exterior Dimensions (W*D*H)(mm)	Lid Opening Height(mm)	Net Weight (kg)	Packing Dimensions (W*D*H)(mm)	Gross Weight(kg)	Certifications
≤58	220~240/50/60	230	270*375*252	526	18	375*560*370	20	CE











BE No.	BE12C5Z00	BE12C4Z00	BE12D4Z00	BE12E4Z00
Туре	Plastic cover	Sealed	Sealed	Sealed
Max Capacity	24*1.5/2ml	24*1.5/2ml	12*5ml	4*8*0.2ml
Max RCF (xg)	25,872	21,382	21,407	18,891
Max Speed (rpm)	16,500	15,000	15,000	15,000
Tubes per Rotor	24	24	12	4
Tube				0000000
Max Tube Capacity (ml)	1.5/2	1.5/2	5	0.2 ml 8-strip PCR
Max Tube Dimensions (Ø x L mm)	11*38	11*38	13.5*75	6*18
Tube Type	Round / Conical bottom	Round /Conical bottom	Round bottom	PCR

<sup>\*</sup>Haier Biomedical also can provide 0.2ml adapter (sets of 24) (BE12J2003), 0.5ml adapter (sets of 24) (BE12JN002) for 1.5/2ml tube rotor





Max. Capacity	24*capillary tubes
Max. RCF (xg)	15,777
Max. Speed (rpm)	12,000



## LX-155T500R

## **Benchtop Multipurpose Refrigerated Centrifuge**

## Scope of Application

It is widely used in fields such as life sciences, biology, medical research, pharmaceuticals, geology, agriculture and forestry, petrochemicals, light industry, and metal particles, for rapid separation and extraction of trace samples.

It also serves as IVD lab equipment for separating human blood, used in diagnostic tests like hematological, immunological screening, and cardiovascular assessment.



## **Product Features**



Short pre-cooling time and precise temperature control performance



Optimized sound design for low noise levels (≤55dB(A))



Wide range of rotors for application flexibility, compatible with tubes/bottles from 0.2 ml to 500ml



Automatic radio-frequency rotor identification (RFID) with over speed protection



Use of brushless induction motors to provide maintenance-free operation, eliminating carbon particle emissions



Real time loading imbalance detection, immediate shutdown protection in case of imbalance



ECO shut-off function reduce energy consumption and extend compressor life



Each rotor use frequency monitoring to guarantee maximum operational safety



**Haier** Biomedical

Intelligent Protection of Life Science

## **Intuitive Interface**



9 acceleration + 10 deceleration



Equipped with user login and password protection, ensuring individualized program data management for each user



Quick recall feature for over 1000+ programs, minimizing time between runs and facilitating quick starts for routine users



The speed, centrifugal force, time and temperature can be changed by user during the centrifugation process



The 7-inch 45° glove-friendly touchscreen, the wide viewing angle more comfortable and practical for users of all heights



## **Product Parts**



Intelligent lid locking monitoring and automatic self-locking to ensure safe operation



Gas springs with safety protection



Drainage tray to prevent internal water accumulation and improve corrosion resistance



Triple-layer construction to ensure safe use



316 stainless steel centrifuge chamber is rust-free and easy



Emergency lid release, pry out the white plastic plug on the side plate with a tool and pull the attached release cord at least 7 times to open the lid



USB port for exporting operation record including speed, time and temperature, operation curve with alarm messages etc

## Key Features of Quick-Flip Biocontainment Safety Lids

## Quick-Flip Hinge Design

Our hinged lids do not need to be removed from the rotor. Just a quick flip gives you fast access to your samples

#### Low-Pressure Closing Clip

The simple low-pressure clip offers a secure closure with an audible "click" and visual cues, enhancing safety and replacing multi-turn screw caps or complex lever systems

#### Ambidextrous & Glove-Friendly

Designed for ease of use, our lids can be operated with either hand and are glove-friendly, ensuring seamless operation

## Premium Polymer Construction

Made from Polyetherimide (PEI), our lids provide excellent chemical resistance, transparency, and durability

#### Certified Safety and Reliability

Independently tested and certified per IEC-61010 Section 2-020 Annex AA, ensuring top-level safety and performance









Model	Maximum Capacity	SpeedRange (rpm)	Max RCF(xg)	Speed Accuracy (rpm)	Time Setting Range	Precooling Function	Temperature Setting Range(°C)	Display
LX-155T500R	4*500 ml	500 ~15,500	26,322	±10	1s~99h59m59s (1-second increment)	Yes	-20 ~40, adjustable to 1	7 inch touch screen

Noise Level dB(A)	Power Supply (V/Hz)	Power(W)	Exterior Dimensions (W*D*H)(mm)	Lid Opening Height(mm)	Net Weight (kg)	Packing Dimensions (W*D*H)(mm)	Gross Weight(kg)	Certifications
≤55	220~240/50/60	1200	658*673*358	860	112	741*799*555	133	CE













BE No.	BE12C9000	BE12CE000	BE12C7000	BE12DS000	BE12DT000	BE12DU000
Туре	Sealed	Sealed	Sealed	Sealed	Sealed	Sealed
Max Capacity	24*1.5/2ml	30*1.5/2ml	48*1.5/2ml	12*10ml	12*15ml	6*50ml
Max RCF (xg)	23,100	26,322	22,132	16,071	21,912	17,065
Max Speed (rpm)	15,500	15,500	14,000	12,500	14,000	12,000
Max Tubes Per Rotor	24	30	48	12	12	6
Tube					темпененный [7]	
Max Tube Capacity (ml)	1.5/2	1.5/2	1.5/2	10	15	50
Max Tube Dimensions (Ø x L mm)	11*38	11*38	11*38	16*86	17*121	29*105
Tube Type	Round / Conical bottom	Round /Conical bottom	Round /Conical bottom	Round bottom	Conical bottom	Round bottom

<sup>\*</sup>Haier Biomedical also can provide 0.2ml adapter (sets of 24) (BE12J2003),0.5ml adapter (sets of 24) (BE12JN002) for 1.5/2ml tube rotor











BE No.	BE12DB000	BE12E5000	BE12DF000	BE12DG000	BE12DE000
Туре	Sealed	Sealed	Steel sleeves	Steel sleeves	Steel sleeves
Max Capacity	6*50ml	12*8*0.2ml	42*15ml	12*50(15)ml	10*100ml
Max RCF (xg)	16,421	14,892	4,472	4,640	3,350
Max Speed (rpm)	12,000	12,000	5,000	5,000	4,500
Max Tubes Per Rotor	6	12	42	12pcs(50ml) 12pcs(15ml)	10
Tube	пинимини)		CIMMINION C	C. Jennessenman	
Max Tube Capacity (ml)	50	0.2 ml 8-strip PCR	15	50(15)	100
Max Tube Dimensions (Ø x L mm)	29*118	6*18	17*121	29*118 17*121	41*120
Tube Type	Conical bottom	PCR	Round/Conical bottom	Round/Conical bottom	Round bottom

#### Adapter(Sets of 6) for BE12DB000

BE No.	BE12JC002	
Max Capacity (ml)	15	
Tube	T-USHIOIMOILO	
Max Tube Dimensions (Ø x L mm)	17*121	
Tube Type	Conical bottom	

## Swing-Out Rotor



Swing-Out Rotor with Round Buckets



BE No.	BE12FJ001	
Max. RCF (xg)	3,184	
Max. speed (rpm)	4.000	



Round Bucket (Sets of 4)

Round Bucket& Lid (sets of 4)

BE No.	BE12G0000		
Max Capacity(ml)	500		
Bucket Type	Flat bottom		



Swing-Out Rotor with **Biocontainment Round Buckets** 



BE No.	BE12FJ001
Max. RCF (xg)	3,184
Max. speed (rpm)	4,000



BE No.	BE12G8000	BE12J2004
Max Capacity (ml)	500	NA
Bucket Type	Flat bottom	NA











BE No.	No Adapter Needed	BE12JP001	BE12FG001	BE12FF001	BE12JS001
Max Capacity	500ml	250ml	225ml	100ml	4*50ml
Tubes Per Adapter/Rotor	4	1/4	1/4	1/4	4/16
Tube					пополинана)
Max Tube Dimension (Ø x L mm)	<b>s</b> 80*136	62*123	60*136	38*120	29*118
Tube Type	Flat bottom	Flat bottom	Conical bottom	Round bottom	Conical bottom

<sup>\*</sup>Haier Biomedical also can provide 0.2ml adapter (sets of 24) (BE12J2003),0.5ml adapter (sets of 24) (BE12JN002) for 1.5/2ml tube rotor









BE No.	BE12JR001	BE12JQ001	BE12JD002	BE12JE002
Max Capacity	10*15ml	14*10ml	19*7/5ml	30*1.5/2ml
Tubes Per Adapter/Rotor	10/40	14/56	19/76	30/120
Tube	mannandi)			
Maximum Tube Dimensions (Ø x L mm)	17*121	16*106	13*75/106	11*38
Tube Type	Conical bottom	Large blood collection tube (Φ16mm)	Short/Long blood collection tube	Round/Conical bottom

<sup>\*</sup>Haier Biomedical also can provide 0.2ml adapter (sets of 24) (BE12J2003),0.5ml adapter (sets of 24) (BE12JN002) for 1.5/2ml tube rotor



Swing-Out Rotor

with Plate Buckets

Adapter (Sets of 4) for BE12G0000,BE12G8000

Swing-Out Rotor Block







Max. RCF (xg)	2,612
Max. speed (rpm)	4,000

BE No. Max Capacity 1 plate **Bucket Type** Rectangular bottom

#### Plates for BE12H5000








Plates per rotor	4	4	4
W*D*H (mm)	128*86*15	128*86*17.5	128*86*20
Plate Type	MTP (without lid)	MTP (with lid)	PCR plate 96-fold


## Plate Swing-out Rotor

BE12FQZ01	
2,381	
4,000	
-	
Rectangular bottom Max height 60mm	
	2,381 4,000





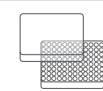

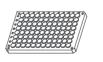





Plates per bucket/rotor	4/8	3/6	3/6	1/2	1/2
W*D*H (mm)	128*86*15	128*86*17.5	128*86*20	128*86*44.5	128*86*60
Plate Type	MTP (without lid)	MTP (with lid)	PCR plate 96-fold	DWP	Filter plate



## **Benchtop High-Speed Refrigerated Centrifuge**



## Scope of Application

It is widely used in fields such as life sciences, biology, medical research, pharmaceuticals, geology, agricul ture and forestry, petrochemicals, light industry, and metal particles, for rapid separation and extraction of trace samples.

It also serves as IVD lab equipment for separating human blood, used in diagnostic tests like hematological, immunological screening, and cardiovascular assessment.

## LX-185T100R





Wide range of rotors for application flexibility, compatible with tubes/bottles from 0.2 ml to 100ml



Short pre-cooling time and precise temperature control performance



Fast sample processing at maximum speeds of up to 18,500 rpm (30,130xg)



Use of brushless induction motors to provide maintenance-free operation, eliminating carbon particle emissions



Optimized sound design for low noise levels (≤52dB(A))



Automatic radio-frequency rotor identification (RFID) with over speed protection



Each rotor use frequency monitoring to guarantee maximum operational safety



Real time loading imbalance detection, immediate shutdown protection in case of imbalance









The 7-inch 45 ° glove-friendly touchscreen, the wide viewing angle more comfortable and practical for users of all heights



Further enhanced with access controls such as user login with password protection, the independent program data of each user is managed separately



Quick recall of up to 1000+ programs with alpha-numeric naming to minimize time between runs and for quick run start by routine users



The speed, centrifugal force, time and temperature can be changed by user during the centrifugation process



9 acceleration + 10 deceleration stages



## Product Parts



Intelligent lid locking monitoring and automatic self-locking to ensure safe operation



Gas springs with safety protection





Triple-layer construction to ensure safe use



316 stainless steel centrifuge chamber is rust-free and easy to clean



Emergency lid release, pry out the white plastic plug on the side plate with a tool and pull the attached release cord at least 7 times to open the lid



Drainage tray to prevent internal water accumulation and improve corrosion resistance



USB port for exporting operation record including speed, time and temperature, operation curve with alarm messages etc



Model	Maximum Capacity	Speed Range (rpm)	Max RCF(xg)	Speed Accuracy (rpm)		Precooling Function	Temperature Setting Range(°C)	Display
LX-185T100R	6*100 ml	500~18500	30,130	±10	1s~99h59m59s (1-second increment)	Yes	-20 ~ 40, adjustable to 1	7 inch touch screen

Noise Level dB(A)	Power Supply (V/Hz)	Power(W)	Exterior Dimensions (W*D*H)(mm)	Lid Opening Height(mm)	Net Weight (kg)	Packing Dimensions (W*D*H)(mm)	Gross Weight(kg)	Certifications
≤52	220~240/50/60	1200	420*680*330	860	57	508*806*528	73	CE





BE No.	BE12DKZ00	BE12C6000	BE12CB000	BE12DJZ00	BE12CC000	BE12CA000
Туре	Sealed	Sealed	Sealed	Sealed	Sealed	Sealed
Max Capacity	12*1.5/2ml	24*1.5/2ml	48*1.5/2ml	12*5ml	12*10ml	8*15ml
Max RCF (xg)	25,636	30,130	19,083	20,124	15,547	14,972
Max Speed (rpm)	18,500	17,500	13,000	15,000	12,500	12,000
Tubes per Rotor	12	24	48	12	12	8
Tube						To the second se
Max Tube Capacity (ml)	1.5/2ml	1.5/2ml	1.5/2ml	5ml	10ml	15ml
Max Tube Dimensions (Ø x L mm)	11*38	11*38	11*38	13.5*75	16*86	17*121
Tube Type	Round /Conical bottom	Round / Conical bottom	Round /Conical bottom	Round bottom	Round bottom	Conical bottom

<sup>\*</sup>Haier Biomedical also can provide 0.2ml adapter (sets of 24) (BE12J2003),0.5ml adapter (sets of 24) (BE12JN002) for 1.5/2ml tube rotor



BE No.	BE12CD000	BE12DD000	BE12DLZ00	BE12DPZ00	BE12DMZ00	BE12DNZ00
Туре	Sealed	Sealed	Sealed	Steel sleeve	Steel sleeve	Steel sleeve
Max Capacity	6*50ml	6*50ml	12*8*0.2ml	18*15ml	6*50/15ml	6*100ml
Max RCF (xg)	15,616	16,421	14,891	3,577	3,577	3,348
Max Speed (rpm)	12,000	12,000	12,000	5,000	5,000	4,800
Tubes per Rotor	6	6	12	18	6 pcs(50ml) 6 pcs(15ml)	6
Tube		шининанф)		THE PROPERTY OF THE PROPERTY O		
Max Tube Capacity (ml)	50	50	0.2ml 8-strip PCR	15	50(15)	100
Max Tube Dimensions (Ø x L mm)	29*105	29*118	6*18	17*121	29*118 17*121	41*120
Tube Type	Round bottom	Conical bottom	PCR	Round /Conical bottom	Round /Conical bottom	Round bottom

#### Adapter(Sets of 6) for BE12DD000

BE No.	BE12JC002	
Max Capacity (ml)	15	
Tube	Дениченный I.	
Max Tube Dimensions (Ø x L mm)	17*121	
Tube Type	Conical bottom	







Swing-Out Rotor Block



BE No.	BE12FD001		
Max. RCF (xg)	3,347		
Max. speed (rpm)	4,600		



Bucket (Sets of 4)



BE No.	BE12H1000
Max Capacity	5*10ml
Tubes per bucket/rotor	5/20
Tube	
Max Tube Dimensions (Ø x L mm)	16*106
Tube Type	Large blood collection tube (Φ16mm)





4,600

BE No.	BE12FD001
Max. RCF (xg)	3,347





BE No.	BE12H2000
Max Capacity	2*15ml
Tubes per bucket/rotor	2/8
Tube	To the state of th
Max Tube Dimensions (Ø x L mm)	17*121
Tube Type	Conical bottom

#### Adapter (Sets of 50) for BE12H1000

BE No.	BE12JN001	
Max Capacity (ml)	5	
Tube		
Max Tube Dimensions (Ø x L mm)	13*75	
Tube Type	Short blood collection tube	U





Max. speed (rpm)



4,600





BE No.	BE12H3000
Max Capacity (ml)	50
Tubes per bucket/rotor	1/4
Tube	minuted)
Max Tube Dimensions (Ø x L mm)	29*118
Tube Type	Conical bottom













NA

NA

BE No.	BE12FD001	4	BE No.	BE12G7000	ВЕ
Max. RCF (xg)	3,347		Max Capacity (ml	100	
Max. speed (rpm)	4,600		Bucket Type	Flat bottom	


















BE No.	BE12JU003	BE12JV003	BE12JW003*	BE12JX003	BE12JY003	BE12JZ003	BE12J0004
Max Capacity	100ml	50ml	2*15ml	15ml	4*10ml	4*5/7ml	4*1.5/2ml
Tubes per adapter/rotor	1/4	1/4	2/8	1/4	4/16	4/16	4/16
Tube		(податилний)	Ct.)Introduction	The second secon			
Max Tube Dimensions (Ø x L mm)	41*120	29*118	17*121	17*121	16*106	13*75/106	11*38
Tube Type	Round bottom	Conical bottom	Conical bottom	Conical bottom	Large blood collection tube (Φ16mm)	Short/Long blood collection tube	Round/Conica bottom

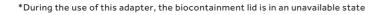





BE12FD001
3,347
4,600



BD CPT Tube Bucket (Sets of 4)



BE No.	BE12G6000
Max Capacity (ml)	8
Tubes per bucket/rotor	1/4
Tube	
Max Tube Dimensions (Ø x L mm)	16*130
Tube Type	Round bottom


BE No.	BE12FC001	
Max. RCF (xg)	2,245	
Max. speed (rpm)	4,500	
Max Capacity	2 plates	
Bucket Type	Rectangular bottom	



Plate Swing-Out Rotor





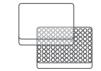





Plates per rotor	4	4	4	2
W*D*H (mm)	128*86*15	128*86*17.5	128*86*20	128*86*44.5
Plate Type	MTP (without lid)	MTP (with lid)	PCR plate 96-fold	DWP





# LX-185T100-J

# **Benchtop High-Speed Centrifuge**

## Scope of Application

Mainly applied to the experimental processes of collecting, separating, and purifying samples such as cells, macromolecular proteins, body fluids, DNA/RNA, plasmids, and viruses. Widely used in universities, third party testing organizations, disease prevention, medical and health, biological pharmaceutical and other related fields.



### **Product Features**



Wide range of rotors for application flexibility, compatible with tubes/bottles from 0.2 ml to 100ml.



Optimized sound design for low noise levels ( $\leq$  63dB(A)).



Fast sample processing at maximum speeds of up to 18,500 rpm (30,130 xg).



Automatic rotor identification locks the maximum rotational speed to prevent overspeed operation.



Use of brushless induction motors to provide maintenance-free operation, eliminating carbon particle emissions.



Real time loading imbalance detection, immediate shutdown protection in case of imbalance.



#### Intuitive Interface



LCD keypad and control knob for easy operation, the wide viewing angle more comfortable and practical for users of all heights.



9 acceleration + 10 deceleration stages.



Quick recall of up to 99 programs for quick run start by routine users.



The speed, centrifugal force, and time can be changed during the centrifugation process

09

Haier Biomedical

Intelligent Protection of Life Science

Benchtop High-Speed Centrifuge

LX-185T100-J

Target speed/centrifugal force/time and real-time operation data display.

Screen lock function to prevent press errors.



One key to recall stored programs for more convenient operation.

Acceleration/Decelera tion rates.



Patented anti-seismic door lock. In case of abnormal vibration, the cover can be locked reliably.



Unique air duct structure design, the heat inside the centrifugal chamber can dissipate quickly.



Dual gas springs for easy door opening and with safety protection.



Triple-layer construction to ensure safe use.



316 stainless steel centrifuge chamber is rust-free and easy to clean.



Emergency lid release, pry out the white plastic plug on the side plate with a tool and pull the attached release cord to open the lid.



Model	Maximum Capacity	Speed Range (rpm)	Max RCF (xg)	Speed Accuracy (rpm)	Time Setting Range	Display	Noise dB (A)	Power Supply (V/Hz)
LX-185T100-J	6*100ml	500-18,500	30,130	±10	1s-99h/59min (1s/1minincrement)	LCD	≤63	220~240/50/60

Power	W)	Exterior Dimensions (W*D*H) (mm)	Lid Opening Height (mm)	Net Weight (kg)	Packing Dimensions (W*D*H)(mm)	Gross Weight(kg)	Certification
500		440*600*346	795	47	495*670*380	51	CE







BE No.	BE12CK000	BE12CL000	BE12CH000	BE12DF001	BE12DH001	BE12DJ001
Туре	Sealed	Sealed	Sealed	Sealed	Sealed	Sealed
Max Capacity	12*1.5/2ml	24*1.5/2ml	48*1.5/2ml	12*5ml	12*10ml	8*15ml
Max RCF (xg)	25,636	30,130	19,083	20,124	15,547	14,972
Max Speed (rpm)	18,500	17,500	13,000	15,000	12,500	12,000
Tubes per Rotor	12	24	48	12	12	8
Tube						To the second se
Max Tube Capacity (ml)	1.5/2ml	1.5/2ml	1.5/2ml	5ml	10ml	15ml
Max Tube Dimensions (Ø xL mm)	11*38	11*38	11*38	13.5*75	16*86	17*121
Tube Type	Round / Conical bottom	Round / Conical bottom	Round / Conical bottom	Round bottom	Round bottom	Conical botton













BE No.	BE12DG001	BE12DK001	BE12E6000	BE12DE001	BE12DD001	BE12DC001
Туре	Sealed	Sealed	Sealed	Steel sleeve	Steel sleeve	Steel sleeve
Max Capacity	6*50ml	6*50ml	12*8*0.2ml	18*15ml	6*50/15ml	6*100ml
Max RCF (xg)	15,616	16,421	14,891	3,577	3,577	3,348
Max Speed (rpm)	12,000	12,000	12,000	5,000	5,000	4,800
Tubes per Rotor	6	6	12	18	6 pcs(50ml) 6 pcs(15ml)	6
Tube		шининий		Tonnen (T)		
Max Tube Capacity (ml)	50ml	50ml	0.2ml 8-strip PCR	15ml	50(15)ml	100ml
Max Tube Dimensions (Ø xL mm)	29*105	29*118	6*18	17*121	29*118 17*121	41*120
Tube Type	Round bottom	Conical bottom	PCR	Round / Conical bottom	Round / Conical bottom	Round bottom

#### Adapter(Sets of 6) for BE12DK001

BE No.	BE12JC002	
Max Capacity (ml)	15	
Tube	Community of the state of the s	
Max Tube Dimensions (Ø x L mm)	17*121	
Tube Type	Conical bottom	









Bucket (Sets of 4)

BE No.	BE12H1000
Max Capacity	5*10ml
Tubes per bucket/rotor	5/20
Tube	
Max Tube Dimensions (Ø x L mm)	16*106
Tube Type	Large blood collection tube (Φ16mm

#### Adapter (Sets of 50) for BE12H1000

BE No.	BE12JN001	
Max Capacity (ml)	5	
Tube	Ī	
Max Tube Dimensions (Ø x L mm)	13*75	
Tube Type	Small blood collection tube	- U





BE No.	BE12FR001
Max. RCF (xg)	3,347
Max. speed (rpm)	4.600





BE No.	BE12H2000
Max Capacity	2*15ml
Tubes per bucket/rotor	2/8
Tube	To the state of th
Max Tube Dimensions (Ø x L mm)	17*121
Tube Type	Conical bottom







Swing-Out Rotor Block



BE No.	BE12FR001
Max. RCF (xg)	3,347
Max. speed (rpm)	4,600





BE No.	BE12H3000
Max Capacity (ml)	50
Tubes per bucket/rotor	1/4
Tube	писаной)
Max Tube Dimensions (Ø x L mm)	29*118
Tube Type	Conical bottom



Swing-Out Rotor Block



BE No.	BE12FR001
Max. RCF (xg)	3,347
Max. speed (rpm)	4,600



BE No.	BE12G7000	BE12J1004
Max Capacity (m	<b>100</b>	NA
Bucket Type	Flat bottom	NA











Bucket (Sets of 4)





BE No.	BE12JU003	BE12JV003	BE12JW003*	BE12JX003	BE12JY003	BE12JZ003	BE12J0004
Max Capacity	100ml	50ml	2*15ml	15ml	4*10ml	4*5/7ml	4*1.5/2ml
Tubes per adapter/rotor	1/4	1/4	2/8	1/4	4/16	4/16	4/16
Tube		ntoround)	CI Jessessesses				
Max Tube Dimensi (Ø x L mm)	ions 41*120	29*118	17*121	17*121	16*106	13*75/106	11*38
Tube Type	Round bottom	Conical bottom	Conical bottom	Conical bottom	Large blood collection tube (Φ16mm)	Short/long blood collection tube	Round/Conical bottom

\*During the use of this adapter, the biocontainment lid is in an unavailable state.

# Swing-Out Rotor





BE No.	BE12FR001	
Max. RCF (xg)	3,347	
Max. speed (rpm)	4,600	






BE No.	BE12G6000
Max Capacity (ml)	8
Tubes per bucket/rotor	1/4
Tube	
Max Tube Dimensions (Ø x L mm)	16*130
Tube Type	Round bottom

BE No.	BE12CJ000	
Max. RCF (xg)	2,245	10 1
Max. speed (rpm)	4,500	60000
Max Capacity	2 plates	
Bucket Type	Rectangular bottom	Plate Swing-Out Rotor



Plates for

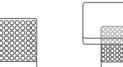





Plates per rotor	4	4	4	2
W*D*H (mm)	128*86*15	128*86*17.5	128*86*20	128*86*44.5
Plate Type	MTP (without lid)	MTP (with lid)	PCR plate 96-fold	DWP



# Benchtop Low-Speed Centrifuge Scope of Application It is widely used in fields such as life sciences, biology, medical research, pharmaceuticals, geology, agriculture

LX-60T100-J

and forestry, petrochemicals, light industry, and metal particles, for rapid separation and extraction of trace samples.

It also serves as IVD lab equipment for separating human blood, used in diagnostic tests like hematological, immunological screening, and cardiovascular assessment.





Offer a compact footprint and compatible with 5ml, 10ml, 15ml, 50ml or 100ml sample tube, small and large blood collection tubes



Use of brushless induction motors to provide maintenance-free operation, eliminating carbon particle emissions



Optimized sound design for low noise levels (≤56dB(A))



Automatic rotor identification with over speed protection



Real time loading imbalance detection, immediate shutdown protection in case of imbalance



10
Benchtop Low-Speed
Centrifuge

LX-60T100-J

Haier Biomedical





LCD keypad and control knob for easy operation, the wide viewing angle more comfortable and practical for users of all heights



Quick recall of up to 99 programs for quick run start by routine users



9 acceleration + 10 deceleration stages



The speed, centrifugal force, and time can be changed during the centrifugation process

Target speed/centrifugal force/time and real-time operation data display



Screen lock function to prevent press errors One key to recall stored programs for more convenient operation

Acceleration/Deceleration rates





Patented anti-seismic door lock. In case of abnormal vibration, the cover can be locked reliably



Gas springs with safety protection





Unique air duct structure design, the heat inside the centrifugal chamber can dissipate quickly



316 stainless steel centrifuge chamber is rust-free and easy to clean



Triple-layer construction to ensure safe use



Emergency lid release, pry out the white plastic plug on the side plate with a tool and pull the attached release cord to open the lid





Model	Maximum Capacity	Speed Range (rpm)	Max RCF(xg)	Speed Accuracy (rpm)	Time Setting Range	Display	Noise Level dB(A)	Power Supply (V/Hz)
LX-60T100-J	6*100ml	500~6,000	4749	±10	1s-99h/59min (1s/1min increment)	LCD	≤56	220~240/50/60

Power(W)	Exterior Dimensions (W*D*H)(mm)	Lid Opening Height (mm)	Net Weight (kg)	Packing Dimensions (W*D*H)(mm)	Gross Weight (kg)	Certifications
175	340*447*281	627	25	441*544*421	28	CE















BE No.	BE11D6000	BE11D5000	BE11DT000	BE11D7000	BE12D6000	BE11DU000
Туре	Steel sleeve	Steel sleeve	Steel sleeve	Steel sleeve	Steel sleeve	Steel sleeve
Max Capacity	6*50(15)ml	12*15ml 12*10ml	18*15ml	6*100ml	12*10ml	36*10ml
Max RCF (xg)	3,578	4,749	3,578	3,349	2,683	2,465
Max Speed (rpm)	5,000	6,000	5,000	4,800	5,000	4,200
Max Tubes Per Rotor	6 pcs(50ml) 6 pcs(15ml)	12	18	6	12	36
Tube			200000000000000000000000000000000000000			
Max Tube Capacity (ml)	50&15ml	15 10ml blood tube	15	100	10	10
Max Tube Dimensions (Ø x L mm)	29*118 17*121	17*121 16*106	17*121	41*120	15*106	15*106
Tube Type	Round/Conical bottom	Round/Conical bottom Large blood collection tube ( <b>Φ</b> 16mm)	Round/Conical bottom	Round bottom	Large blood collection tube ( <b>Φ</b> 15mm)	Large blood collection tube (Φ15mm

#### Adapter(Sets of 50) for BE12D6000, BE11DU000

BE12JN001	
5	
13*75	
Small blood collection tube	W
	13*75









Swing-Out Rotor Block





BE No.	BE12H8000	
Max Capacity	1*15ml 2*10ml 2*5ml	
Tubes per bucket/rotor	15ml(1/4),10ml(2/8) 5ml(2/8)	
Tube		
Max Tube Dimensions (Ø x L mm)	17*121 16*106 13*75	
Tube Type	15ml(Conical bottom) 10ml/5ml(Round bottom)	







BE No.	BE12F2001
Max. RCF (xg)	2,770
Max. speed (rpm)	4,200



t of 4)	
· .,	

BE No.	BE12H1000
Max Capacity	5*10ml
Tubes per bucket/rotor	5/20
Tube	
Max Tube Dimensions (Ø x L mm)	16*106
Tube Type	Large blood collection tube ( <b>Φ</b> 16mm)







BE No.	BE12H3000
Max Capacity (ml)	50
Tubes per bucket/rotor	1/4
Tube	(пешено)
Max Tube Dimensions (Ø x L mm)	29*118
Tube Type	Conical bottom

Adapter	Sets of	50) for B	E12H1000
Auaptei	JELS UI .	ם וטו נטכ	LIZITIO

BE No.	BE12JN001	
Max Capacity (ml)	5	
Tube		
Max Tube Dimensions (Ø x L mm)	13*75	
Tube Type	Small blood collection tube	V

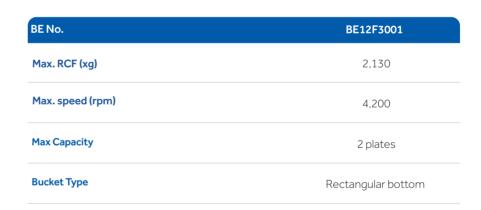




Plate Swing-out Rotor





BE No.	BE12F2001
Max. RCF (xg)	2,770
Max. speed (rpm)	4,200





BE No.	BE12H2000
Max Capacity	2*15ml
Tubes per bucket/rotor	2/8
Tube	- Пенення (D
Max Tube Dimensions (Ø x L mm)	17*121
Tube Type	Conical bottom











Plates per rotor	4	4	4	2
W*D*H (mm)	128*86*15	128*86*17.5	128*86*20	128*86*44.5
Plate Type	MTP (without lid)	MTP (with lid)	PCR plate 96-fold	DWP



# Benchtop Low-Speed Large-Capacity Centrifuge



LX-60T500-J

# Scope of Application

It is widely used in fields such as life sciences, biology, medical research, pharmaceuticals, geology, agriculture and forestry, petrochemicals, light industry, and metal particles, for rapid separation and extraction of trace samples.

It also serves as IVD lab equipment for separating human blood, used in diagnostic tests like hematological, immunological screening, and cardiovascular assessment.





Wide range of rotors for application flexibility, compatible with tubes/bottles from 0.2 ml to 500ml



Automatic rotor identification with over speed protection



Use of brushless induction motors to provide maintenance-free operation, eliminating carbon particle emissions



Real time loading imbalance detection, immediate shutdown protection in case of imbalance



Optimized sound design for low noise levels (≤63dB(A))





LCD keypad and control knob for easy operation, the wide viewing angle more comfortable and practical for users of all heights



Quick recall of up to 99 programs for quick run start by routine users



9 acceleration + 10 deceleration stages



The speed, centrifugal force, and time can be changed during the centrifugation process





Product Parts

Patented anti-seismic door lock. In case of abnormal vibration, the cover can be locked reliably





Unique air duct structure design, the heat inside the centrifugal chamber can dissipate quickly



Triple-layer construction to ensure safe use





Emergency lid release, pry out the white plastic plug on the side plate with a tool and pull the attached release cord to open the lid



316 stainless steel centrifuge chamber is rust-free and easy to clean



## Key Features of Quick-Flip Biocontainment Safety Lids

#### Quick-Flip Hinge Design

Our hinged lids do not need to be removed from the rotor. Just a quick flip gives you fast access to your samples

#### Low-Pressure Closing Clip

The simple low-pressure clip offers a secure closure with an audible "click" and visual cues, enhancing safety and replacing multi-turn screw caps or complex lever systems

#### Ambidextrous & Glove-Friendly

Designed for ease of use, our lids can be operated with either hand and are glove-friendly, ensuring seamless operation

#### Premium Polymer Construction

Made from Polyetherimide (PEI), our lids provide excellent chemical resistance, transparency, and durability

#### Certified Safety and Reliability

Independently tested and certified per IEC-61010 Section 2-020 Annex AA, ensuring top-level safety and performance

Choose our Quick-Flip biocontainment safety lids for innovative, safe, and durable lab solutions.



Model	Maximum Capacity	Speed Range (rpm)	Max RCF(xg)	Speed Accuracy (rpm)	Time Setting Range	Display	Noise Level dB(A)	Power Supply
LX-60T500-J	4*500 ml	500~6000	5353	±10	1s-99h/59min (1s/1min increment)	LCD	≤63	220~240/50/60

Power(W)	Exterior Dimensions (W*D*H)(mm)	Lid Opening Height (mm)	Net Weight (kg)	Packing Dimensions (W*D*H)(mm)	Gross Weight (kg)	Certifications
700	478*624*365	840	60	560*752*562	75	CE
















BE No.	BE12DA000	BE12D7000	BE12C8000	BE12D8000	BE12D9000	BE12DQ000
Туре	Steel sleeve	Steel sleeve	Steel sleeve	Steel sleeve	Steel sleeve	Steel sleeve
Max Capacity	6*50 (15)ml	12*50 (15)ml	12*15ml	42*15ml	10*100ml	60*10ml
Max RCF (xg)	5,353	4,528	5,353	4,528	3,260	2,576
Max Speed (rpm)	6,000	5,000	6,000	5,000	4,500	4,000
Max Tubes Per Rotor	6 pcs(50ml) 6 pcs(15ml)	12 pcs(50ml) 12 pcs(15ml)	12	42	10	60
Tube			Processor (1)			
Max Tube Capacity (ml)	50 (15)ml	50 (15)ml	15ml	15ml	100ml	10ml
Max Tube Dimensions (Ø x L mm)	29*118 17*121	29*118 17*121	17*121	17*121	41*120	15*106
Tube Type	Round/Conical bottom	Round/Conical bottom	Round/Conical bottom	Round/Conical bottom	Round bottom	Large blood collection tube ( \$\Pi\$ 15mm)

#### Adapter(Sets of 50) for BE12DQ000

BE No.	BE12JN001	
Max Capacity (ml)	5	
Tube		
Max Tube Dimensions (Ø x L mm)	13*75	
Tube Type	Small blood collection tube	







with round buckets





Round Bucket Sets of 4)	
----------------------------	--



	BE No.	BE12F
	Max. RCF (xg)	3,1
Swing-out rotor	May speed (rpm)	4.0

Swing-Out Rotor Block

BE No.	BE12FK001	
Max. RCF (xg)	3,184	_1
Max. speed (rpm)	4,000	7

BE No.	BE12G0000
Max Capacity (ml)	500
Bucket Type	Flat bottom



Swing-Out Rotor with Biocontainment Round Buckets



BE No.	BE12FK001	
Max. RCF (xg)	3,184	
Max. speed (rpm)	4,000	





BE No.	BE12G8000	BE12J2004
Max Capacity (ml)	500	NA
Bucket Type	Flat bottom	NA









BE No.	No Adapter Needed	BE12JP001	BE12FG001	BE12FF001	BE12JS001
Max Capacity	500ml	250ml	225ml	100ml	4*50ml
Tubes Per Adapter/Rotor	4	1/4	1/4	1/4	4/16
Tube					полинице))
Max Tube Dimensions (Ø x L mm)	80*136	62*123	60*136	38*120	29*118
Tube Type	Flat bottom	Flat bottom	Conical bottom	Round bottom	Conical bottom











BE No.	BE12JR001	BE12JQ001	BE12JD002	BE12JE002
Max Capacity	10*15ml	14*10ml	19*5/7ml	30*1.5/2ml
Tubes Per Adapter/Rotor	10/40	14/56	19/76	30/120
Tube	(houseoned			
Maximum Tube Dimensions (Ø x L mm)	17*121	16*106	13*75/106	11*38
Tube Type	Conical bottom	Large blood collection tube( <b>Φ</b> 16mm)	Short/Long blood collection tube	Round/Conical bottom

<sup>\*</sup>Haier Biomedical also can provide 0.2ml adapter (sets of 24) (BE12J2003), 0.5ml adapter (sets of 24) (BE12JN002) for 1.5/2ml tube rotor.





Swing-Out Rotor with Conical Buckets







BE No.	BE12FK001	
Max. RCF (xg)	3,506	
Max. speed (rpm)	4,000	

BE No.	BE12G1000	
Max Capacity (ml)	500	
Bucket Type	Conical bottom	

# Adapter (Sets of 4) for BE12G1000



BE No.	No adapter needed	BE12JB002
Max Capacity (ml)	500	250
Tubes Per Adapter/Rotor	4	1/4
Tube		
Max Tube Dimensions (Ø x L mm)	96*150	60*165
Tube Type	Conical bottom	Conical bottom



Swing-out Rotor Block

Max. RCF (xg)

Max. speed (rpm)

BE No.

Swing-Out Rotor with Square Buckets



BE12FK001	
3,184	
4 000	



	BE No.	BE12HA000			
	Max Capacity (ml)	500			
	Bucket Type	Flat bottom			

Square Bucket (Sets of 4)

# Decap Adapter (Sets of 4) for BE12HA000



BE No.	BE12JV001	
Max Capacity	24*5/7ml	
Tubes Per Adapter/Rotor	24/96	
Tube		
Max Tube Dimensions (Ø x L mm)	13*75/106	
Tube Type	Short/Long blood collection tube	







BE No.	BE12FK001			
Max. RCF (xg)	2,612			
Max. speed (rpm)	4,000			



**Bucket Type** 











Rectangular bottom

Plates per rotor	4	4	4
W*D*H (mm)	128*86*15	128*86*17.5	128*86*20
Plate Type	MTP (without lid)	MTP (with lid)	PCR plate 96-fold

# **Mini-Centrifuge**



**LX-120T2Z / 50T2Z** 

# Scope of Application

This product is used in micro centrifugal scenarios (single tube capacity is less than or equal to 2ml)











LX-50T2Z LX-120T2Z



Haier Biomedical Intelligent Protection of Life Science





#### Quiet

Suspension with shock absorber reduces noise output



#### Convenient

Simple and fast operation



#### Small

Can be held in your hand



#### Safe

- Safety interlock lid with automatic shut-off if lid is open
- Anti vibration & absorption suction cup
- Splash and leak protection



Model	LX-120T2Z	LX-50T2Z
Maximum Capacity	8*1.5/2 ml	8*1.5/2 ml
Max Speed (rpm)	12,000	5,000
Max RCF (xg)	7,325	1,270
Time Setting Range	10s-99m50s	-
Display	OLED	-
Noise (dB(A))	≤50	≤40
Power Supply (V/Hz)	100-240/50/60	100-240/50/60
Power (W)	35	25
Exterior Dimensions (W*D*H)(mm)	152*142*122	152*142*122
Net Weight (kg)	1.2	1.2
Packing Dimensions (W*D*H)(mm)	265*205*200	265*205*200
Gross Weight (kg)	3	3







Туре	Fixed-Angle Rotor	Fixed-Angle Rotor
Max Capacity	8*1.5/2 ml	2*8*0.2 ml
Max RCF (xg)	7,325	6,225
Max Speed (rpm)	12,000	12,000
Tubes per Rotor	8	2
- Tube		9000000
Max Tube Capacity(ml)	1.5/2	0.2 8-strip PCR
Max Tube Dimension (Ø x L mm)	11*38	6*18
Tube Type	Round/Conical bottom	PCR







Туре	Fixed-Angle Rotor	Fixed-Angle Rotor		
Max Capacity	8*1.5/2 ml	2*8*0.2 ml		
Max RCF (xg)	1,270	1,150		
Max Speed (rpm)	5,000	5,000		
Tubes per Rotor	8	2		
Tube		900000		
Max Tube Capacity (ml)	1.5/2	0.2 8-strip PCR		
Max Tube Dimension (Ø x L mm)	11*38	6*18		
Tube Type	Round/Conical bottom	PCR		

#### Adapter(Sets of 24) for 1.5/2ml Fixed-angle rotor

BE No.	BE12J2003	BE12JN002
Max Capacity (ml)	0.2	0.5
Tube	Ð	
Max Tube Dimensions (Ø x L mm)	6.3*23	8*35
Tube Type	Conical bottom	Conical bottom

Centrifuge Consumables





- Cap: HDPE, Tube: PP
- Temperature Range: -80°C -121°C







CAT.NO.	Volume (mL)	Color	RCF	Sterile	Packing
311005	15	Blue Cap	16000xg	√	25 pcs/bag, 20 bags/ case
311006	50	Blue Cap	16000xg	√	25 pcs/bag, 20 bags/ case















- Holding samples, mixing reagents. Low temperature storage of antibodies, nucleic acids, etc. Sample preparation during the experiment. Centrifuge separation and precipitation of trace samples.
- Temperature Range: -80°C -121°C.

CAT.NO.	Volume (mL)	Material	RCF (xg)	Color	Sterile	Packing
311001	0.6	PP	20000	Clear	×	1000 pcs/box, 10 boxes/ case
311002	1.5	PP	20000	Clear	×	500 pcs/box, 10 boxes/ case
311004	2.0	PP	20000	Clear	×	500 pcs/box, 10 boxes/ case








- Volume Selection: 225/250/500 mL
- Temperature Range: -80°C -121°C


CAT.NO.	Volume (mL)	Туре	Material	RCF (xg)	Sterile	Packing
321007	225	Conical bottom	Tube: PP;Cap: HDPE	6000	$\checkmark$	6 pcs/bag, 6 bags/ case
321008	250	Conical bottom	Tube: PP;Cap: HDPE	6000	√	6 pcs/bag, 17 bags/ case
321009	500	Conical bottom	Tube: PP;Cap: HDPE	6000	√	6 pcs/bag, 6 bags/ case
321010	500	Flat bottom	PP	8000	×	6 pcs/bag, 10 bags/ case



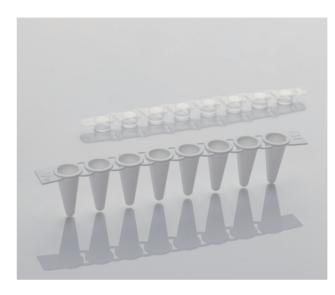



- Volume Selection: 100/200 μL
- Temperature Range: -80°C -121°C

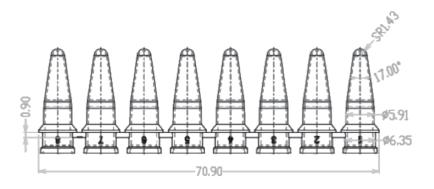

Product	CAT.NO.	Volume (mL)	Color	Product Description	Packing
PCR tubes	211001	0.2	Clear	single tube, thin-wall, flat cap	500pcs/bag, 2 bags/box, 10 boxes/case
	222002	0.1	Clear	8-strip tube, flat cap	125 sets/bag, 10 bags/case
	222003	0.1	White	8-strip tube, flat cap	125 sets/bag, 10 bags/case
	292004	0.1	Clear	8-strip tubetype B, flat cap	125 sets/bag, 10 bags/case
PCR 8-strip	292005	0.1	White	8-strip tubetype B, flat cap	125 sets/bag, 10 bags/case
tubes	292006	0.2	Clear	8-strip tubetype B, flat cap	125 sets/bag, 10 bags/case
	292007	0.2	White	8-strip tubetype B, flat cap	125 sets/bag, 10 bags/case
	222008	0.2	Clear	8-strip tube, flat cap	125 sets/bag, 10 bags/case
	272009	0.2	Clear	8-strip tube, convex cap	125 sets/bag, 10 bags/case







The PCR series products (Polymerase Chain Reaction, PCR) are mainly used as primers, Taq DNA polymerase, dNTP, template nucleic acid, Mg, buffer, and other carriers involved ,They are widely used in genetics, biochemistry, immunity, medicine, and other fields. They are not only used in basic research such as gene isolation, cloning, and nucleic acid sequence analysis, but also in disease diagnosis or anywhere with DNA or RNA.

PCR are used in batch samples to meet different user needs and have wide applicability, making them suitable for different autom ation applications.







CAT.NO.	Volume (mL)	Туре	Color	Sterile	Packing
232024	0.2	No Skirt	Clear	×	10pcs/box, 5boxes/case
242025	0.2	Half Skirt	Clear	×	10pcs/box, 5boxes/case
242026	0.2	Half Skirt,Anti-static	Clear	×	10pcs/box, 5boxes/case
232029	0.1	No Skirt	Clear	×	10pcs/box, 5boxes/case
242030	0.1	Half Skirt	Clear	×	10pcs/box, 5boxes/case
242031	0.1	Half Skirt,Anti-static	Clear	×	10pcs/box, 5boxes/case
252032	0.1	Full Skirt	Clear	×	10pcs/box, 5boxes/case
252035	0.1	Full Skirt,H1	Clear+White	×	10pcs/box, 5boxes/case
252033	0.1	Full Skirt,Anti-static	Clear	×	10pcs/box, 5boxes/case



Product	CAT.NO.	Dimension	Temperature Range	Packing
	202034	140mm*80mm	-20°C~120°C	100pcs/box